Graphical Symbols
Symboles graphiques
Símbolos Gráficos

Attention: see instructions for use.
Attention: voir le mode d'emploi.
Atención: Vea las instrucciones de uso.

Keep dry.
Conserver au sec.
Consérvelos secos.

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hollisterwoundcare and wave logo are trademarks of Hollister Incorporated. Restore is a trademark of Hollister Wound Care, LLC.

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Manufactured for:
Hollister Wound Care, LLC
Libertyville, Illinois 60048
Distributed in Canada by:
Distribué au Canada par:
Hollister Limited
Aurora, Ontario L4G 1G3

Hydrocolloid Dressing with Foam Backing

Hydrocolloid Dressings are sterile, occlusive dressings. The flexible outer layer is impermeable to micro-organisms, urine and feces, and protects wounds from the external environment. The heat-activated, self-adhesive inner layer maintains a moist environment while absorbing excess wound exudate.

The dressings are designed for use on dermal ulcers and partial-thickness wounds with light to moderate amounts of exudate.

Indications for Use

- Venous stasis ulcers
- Superficial wounds
- Pressure ulcers (Stages I and II)
- Arterial ulcers
- Diabetic ulcers
- Surgical incisions
- Traumatic wounds

Contraindications

- Active vasculitis
- Infection
- Pressure ulcers (Stages III and IV)

Instructions for Use

Preparing the Wound Area

Before applying a dressing, the wound should be gently rinsed or irrigated with an appropriate solution. The periwound skin should be cleaned and dried for secure application.

Applying the Dressing




To ensure proper adhesion, the dressing should extend at least one inch beyond the wound edge. Dressings may be overlapped or cut to accommodate the size or location of the wound.

1. Remove the release paper from the dressing.
2. Center the dressing over the wound, being careful to minimize touching the adhesive side. Press the dressing in place. Initial tack may be improved by warming the dressing with your hands prior to application or after dressing is in place.

Removing the Dressing

3. To remove, carefully lift an edge of the dressing while pressing down on the skin adjacent to the edge. Continue this procedure around the wound until all of the edges are free of the skin. Gently lift the dressing off the wound.

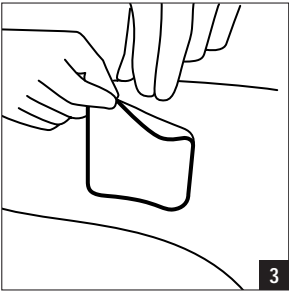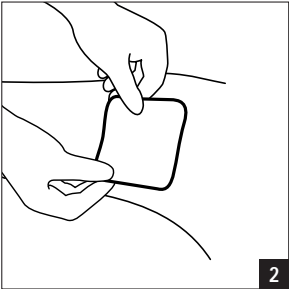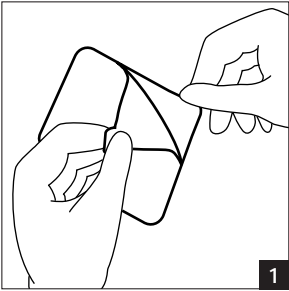Gently rinse or irrigate the wound as needed, remove excess moisture, and apply a new dressing.

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© 2006 Hollister Incorporated.